


# AUTO CHANGE OVER UNIT

Automated Control panel is a panel attached to 2-4 banks having multiple gas cylinders. Automated Control Panel System contains analog pressure meter for displaying pressure of each incoming bank and mainline gas pressure. As soon as the gas in one bank runs below set pressure, it sounds an alarm and switches to the other bank. Specific colored LED lights indicate the bank in use and the bank on standby. The device also indicates whether the gas pressure is in normal or low along the banks and the mains



Semi-Automatic Change-Over Unit

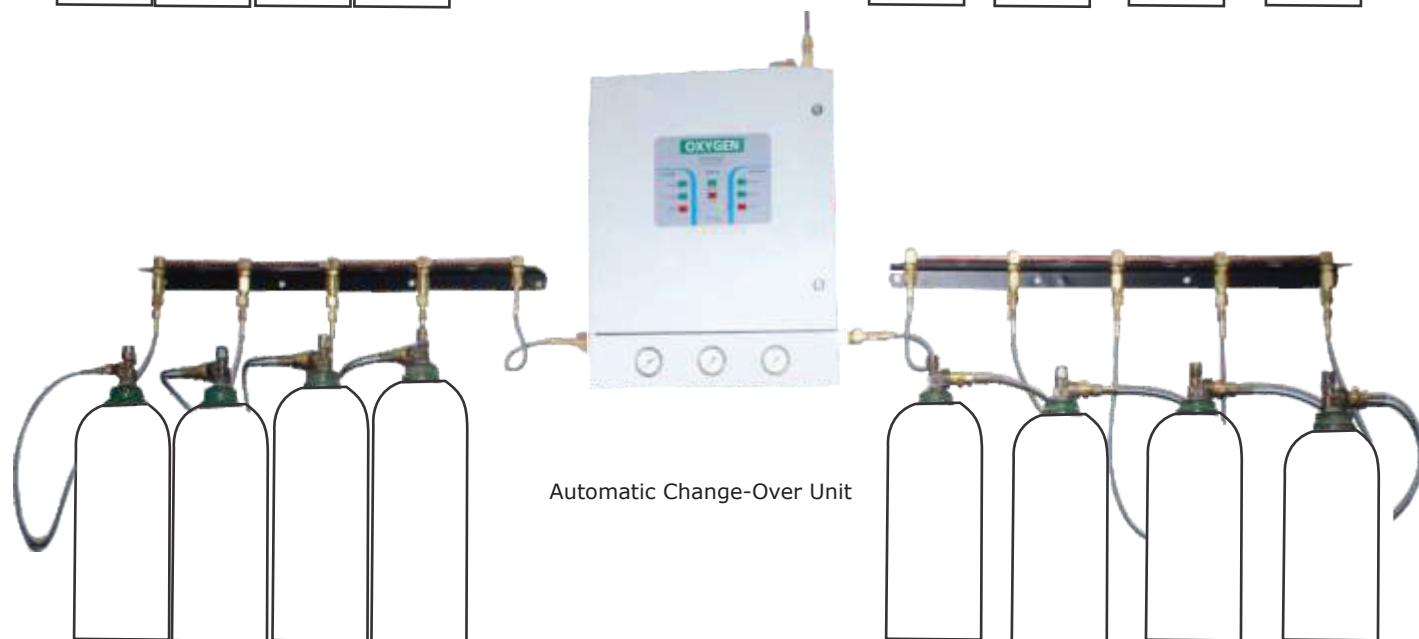

Automatic Change-Over Unit